


SEDEX Low Temperature Evaporative Light Scattering Detector

SAMPLE: AMINO ACIDS

HPLC

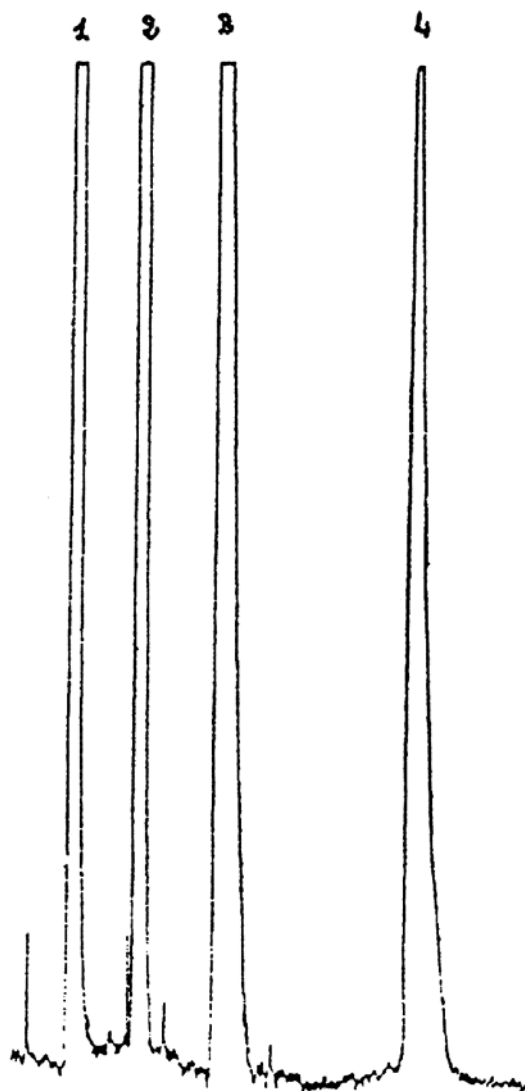

1: Leucine
2: Argenine
3: Leucine
4: Proline

CONDITIONS

Column	<i>Zorbax SB C18</i>
Eluent	Water + 0.3 % Trifluoroacetic Acid
Flow Rate	<i>1.0 ml/min</i>
Pressure	<i>2 bars</i>
Gas	<i>Air</i>
Temperature	<i>40° C</i>

Analyzed by: CE Saclay - Service des Molecules Marquees - FRANCE

S.E.D.E.R.E.
PARC VOLTA - BP 27
9, rue Parmentier
94141 Alfortville Cedex FRANCE

Application #90